Product information

Information about other products is available at: www.demeditec.com

Serotonin ELISA

REF

DEE8900

96

Demeditec Diagnostics GmbH
Lise-Meitner-Strasse 2
24145 Kiel – Germany
www.demeditec.com

CONTENT

| | | |
|-------|---|----|
| 1. | INTRODUCTION | 3 |
| 1.1 | INTENDED USE AND PRINCIPLE OF THE TEST..... | 3 |
| 1.2 | CLINICAL APPLICATION..... | 3 |
| 2. | PROCEDURAL CAUTIONS, GUIDELINES, WARNINGS AND LIMITATIONS..... | 3 |
| 2.1 | PROCEDURAL CAUTIONS, GUIDELINES AND WARNINGS..... | 3 |
| 2.2 | LIMITATIONS | 4 |
| 2.2.1 | INTERFERING SUBSTANCES AND PROPER HANDLING OF SPECIMENS..... | 4 |
| 2.2.2 | DRUG AND FOOD INTERFERENCES..... | 4 |
| 2.2.3 | HIGH-DOSE-HOOK EFFECT | 4 |
| 3. | STORAGE AND STABILITY | 5 |
| 4. | MATERIALS | 5 |
| 4.1 | CONTENTS OF THE KIT | 5 |
| 4.2 | CALIBRATION AND CONTROLS..... | 6 |
| 4.3 | ADDITIONAL MATERIALS REQUIRED BUT NOT PROVIDED IN THE KIT | 6 |
| 4.4 | ADDITIONAL EQUIPMENT REQUIRED BUT NOT PROVIDED IN THE KIT | 6 |
| 5. | SAMPLE COLLECTION AND STORAGE..... | 6 |
| 6. | TEST PROCEDURE..... | 7 |
| 6.1 | PREPARATION OF REAGENTS AND FURTHER NOTES..... | 7 |
| 6.2 | PREPARATION OF SAMPLES – ACYLATION | 7 |
| 6.3 | SEROTONIN ELISA | 8 |
| 7. | CALCULATION OF RESULTS..... | 8 |
| 7.1 | EXPECTED REFERENCE VALUE | 9 |
| 7.2 | TYPICAL STANDARD CURVE | 9 |
| 8. | QUALITY CONTROL | 9 |
| 9. | ASSAY CHARACTERISTICS..... | 10 |
| 9.1 | PERFORMANCE DATA..... | 10 |
| 9.2 | METROLOGICAL TRACEABILITY | 11 |
| 10. | REFERENCES/LITERATURE..... | 12 |
| 11. | CHANGES | 12 |

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1. Introduction

1.1 Intended use and principle of the test

Enzyme Immunoassay for the quantitative determination of serotonin in urine and serum, to evaluate serotonin homeostasis. The determination of serotonin in urine is helpful for the determination of neurostress.

The quantitative determination of serotonin follows the basic principles of the enzyme immunoassay. In the first step, serotonin is quantitatively acylated. The subsequent competitive ELISA kit uses the microtiter plate format. The antigen is bound to the solid phase of the microtiter plate. The acylated standards, controls and samples and the solid phase bound analyte compete for a fixed number of antiserum binding sites. After the system is in equilibrium, free antigen and free antigen-antiserum complexes are removed by washing. The antibody bound to the solid phase is detected by an anti-rabbit IgG-peroxidase conjugate using TMB as a substrate. The reaction is monitored at 450 nm.

Quantification of unknown samples is achieved by comparing their absorbance with a reference curve prepared with known standard concentrations. Manual processing is recommended. The use of automatic laboratory equipment is the responsibility of the user. This in-vitro diagnostic is for professional use only.

1.2 Clinical application

Serotonin (5-hydroxytryptamine) is an intermediate product of tryptophan metabolism [2], a well-studied neurotransmitter, and may also act as a peripheral hormone [3]. Synthesis occurs mainly in enterochromaffin cells (ec-cells) of the gastrointestinal tract and in neurons [2, 4]. It is present in high concentrations in ec-cells of the intestine, serotonergic neurons of the brain, and platelets [1, 2, 4 – 6]. Serotonin is mainly degraded to 5-hydroxyindole acetic acid (5-HIAA) or melatonin [2, 7] and can be excreted in the urine [8]. In the bloodstream, the vast majority of serotonin is found in platelets [9] and can be readily detected in serum. Altered serotonin levels in serum and/or urine can indicate both physical and psychological dysfunction.

Serotonin balance may be impaired in serum and/or urine in a variety of conditions. For example, decreased serotonin levels have been demonstrated in depression, anxiety, and even pain sensitivity compared to unaffected subjects [6, 8, 10]. Increased serotonin levels, on the other hand, have been reported in patients with serotonin-secreting neuroendocrine tumors, also called carcinoid tumors [3, 12, 13], or hepatocellular carcinomas [14].

Therapeutic consequences should never be based on laboratory results alone even if all test results are in agreement with the items as under point "Procedural cautions, guidelines and warnings". Any laboratory result is only a part of the total clinical picture of the patient.

Only in cases where the laboratory results are in an acceptable agreement with the overall clinical picture of the patient it can be used for therapeutic consequences.

The test result itself should never be the sole determinant for deriving any therapeutic consequences.

2. Procedural cautions, guidelines, warnings and limitations

2.1 Procedural cautions, guidelines and warnings

- (1) This kit is intended for professional use only. Users should have a thorough understanding of this protocol for the successful use of this kit. Only the test instruction provided with the kit is valid and must be used to run the assay. Reliable performance will only be attained by strict and careful adherence to the instructions provided.
- (2) This assay was validated for a certain type of sample as indicated in Intended Use (please refer to Chapter 1). Any off-label use of this kit is in the responsibility of the user and the manufacturer cannot be held liable.
- (3) The principles of Good Laboratory Practice (GLP) must be followed.
- (4) To reduce exposure to potentially harmful substances, wear lab coats, disposable protective gloves and protective glasses where necessary.
- (5) If serious incidents should occur in connection with this product, they should be reported to the manufacturer and the competent national authorities.
- (6) All kit reagents and specimens should be brought to room temperature and mixed gently but thoroughly before use. For dilution or reconstitution purposes, use deionized, distilled, or ultrapure water. Avoid repeated freezing and thawing of reagents and specimens.
- (7) The microplate contains snap-off strips. Unused wells must be stored at 2 – 8 °C in the sealed foil pouch with desiccant and used in the frame provided. Microtiter strips which are removed from the frame for usage should be marked accordingly to avoid any mix-up.
- (8) Duplicate determination of sample is highly recommended.
- (9) Once the test has been started, all steps should be completed without interruption. Make sure that the required reagents, materials, and devices are prepared for use at the appropriate time.

-
- (10) Incubation times do influence the results. All wells should be handled in the same order and time intervals.
 - (11) To avoid cross-contamination of reagents, use new disposable pipette tips for dispensing each reagent, sample, standard and control.
 - (12) A standard curve must be established for each run.
 - (13) The controls should be included in each run and fall within established confidence limits. The confidence limits are listed in the QC-Report provided with the kit.
 - (14) Do not mix kit components with different lot numbers within a test and do not use reagents beyond expiry date as shown on the kit labels.
 - (15) Avoid contact with Stop Solution containing 0.25 M H₂SO₄. It may cause skin irritation and burns. In case of contact with eyes or skin, rinse off immediately with water.
 - (16) TMB substrate has an irritant effect on skin and mucosa. In case of possible contact, wash eyes with an abundant volume of water and skin with soap and abundant water.
 - (17) For information about hazardous substances included in the kit please refer to Safety Data Sheet (SDS). The Safety Data Sheet for this product is made available directly on the website of the manufacturer or upon request.
 - (18) Kit reagents must be regarded as hazardous waste and disposed of according to national regulations.
 - (19) The expected reference values reported in this test instruction are only indicative. It is recommended that each laboratory establishes its own reference intervals.
 - (20) In case of any severe damage to the test kit or components, the manufacturer has to be informed in writing, at the latest, one week after receiving the kit. Severely damaged single components must not be used for a test run. They must be stored properly until the manufacturer decides what to do with them. If it is decided that they are no longer suitable for measurements, they must be disposed of in accordance with national regulations.
 - (21) The results obtained with this test kit should not be taken as the sole reason for any therapeutic consequence but must be correlated to other diagnostic tests and clinical observations.

2.2 Limitations

Any inappropriate handling of samples or modification of this test might influence the results.
If you have any further questions, please contact the manufacturer.

2.2.1 Interfering substances and proper handling of specimens

Urine

Please note the sample collection! It cannot be excluded that high acid concentrations lead to incorrect results. Up to 30 µl 100% acetic acid per 1 ml urine no influence on the results was observed.

Serum

Samples containing precipitates or fibrin strands might cause inaccurate results.

Hemolytic samples (up to 2 mg/ml hemoglobin), icteric samples (up to 50 mg/dl bilirubin) and lipemic samples (up to 834 mg/dl triglycerides) have no influence on the assay results.

If the concentrations cannot be estimated and there are doubts as to whether the above limit values for hemolytic, icteric or lipemic samples are complied with, the samples should not be used in the assay.

2.2.2 Drug and food interferences

The following foods and stimulants can affect the serotonin content in the sample. Alcohol, pineapple, eggplant, avocados, bananas, grapefruit, currants, cocoa, kiwis, caffeine, melons, mirabelles, nicotine, pecans, peaches, plums, chocolate, gooseberries, tomatoes, walnuts.

Some drugs can also affect serotonin levels in the sample. For example, taking amphetamines, acetanilide, coumarins, ephedrine, guaifenesin, mephenesin (carbamate), methocarbamol, monoamine oxidase inhibitors (MAO inhibitors), acetaminophen, phenacetin, phenobarbital, phentolamine, or reserpine can lead to increased serotonin levels. In contrast, acetylsalicylic acid, chlorpromazine, isoniazid, levodopa, methenamine, methyldopa, promethazine, selective serotonin reuptake inhibitors (SSRIs), or streptozocin may result in decreased serotonin levels.

Therefore, 2 – 4 days prior to specimen collection, these foods should be avoided and the medications discontinued if medically justifiable.

2.2.3 High-Dose-Hook effect

No hook effect was observed in this test.

3. Storage and stability

Store kit and reagents at 2 – 8 °C until expiration date. Do not use components beyond the expiry date indicated on the kit labels. Once opened, the reagents are stable for 2 months when stored at 2 – 8 °C. Once the resealable pouch has been opened, care should be taken to close it tightly again including the desiccant.

4. Materials

4.1 Contents of the kit

| | | |
|-------------------------|---|---|
| BA E-0030 | WASH-CONC 50x | Wash Buffer Concentrate – concentrated 50x |
| Content: | Buffer with a non-ionic detergent and physiological pH | |
| Volume: | 1 x 20 ml/vial, purple cap | |
| BA E-0040 | CONJUGATE | Antibody Conjugate – ready to use |
| Content: | Goat anti-rabbit immunoglobulins conjugated with peroxidase | |
| Volume: | 1 x 12 ml/vial | |
| Description: | Species is goat, red cap | |
| BA E-0055 | SUBSTRATE | Substrate – ready to use |
| Content: | Chromogenic substrate containing 3,3',5,5'-tetramethylbenzidine, substrate buffer and hydrogen peroxide | |
| Volume: | 1 x 12 ml/vial, black cap | |
| BA E-0080 | STOP-SOLN | Stop Solution – ready to use |
| Content: | 0.25 M sulfuric acid | |
| Volume: | 1 x 12 ml/vial, grey cap | |
| Hazards identification: |  | |
| | | H290 May be corrosive to metals. |
| BA E-0931 | U SER 5-HIAA | Serotonin Microtiter Strips – ready to use |
| Content: | 1 x 96 well (12x8) antibody precoated microwell plate in a resealable white pouch with desiccant | |
| BA E-6612 | ACYL-REAG | Acylation Reagent – ready to use |
| Content: | Acylation reagent in DMSO | |
| Volume: | 2 x 3 ml/vial, white cap | |
| BA E-8910 | SER-AS | Serotonin Antiserum – ready to use |
| Content: | Rabbit anti-Serotonin antibody, blue coloured | |
| Volume: | 1 x 6 ml/vial, blue cap | |
| Description: | Species is rabbit | |
| BA E-8911 | ACYL-BUFF | Acylation Buffer – ready to use |
| Content: | TRIS buffer with non-mercury preservative | |
| Volume: | 1 x 55 ml/vial, grey cap | |

4.2 Calibration and Controls

Standards and Controls – ready to use

| Cat. no. | Component | Colour/Cap | Concentration [ng/ml] (=µg/l) | Concentration [nmol/l] | Volume/Vial |
|-----------|-------------------|------------|---|------------------------|-------------|
| BA R-8901 | STANDARD A | white | 0 | 0 | 4 ml |
| BA R-8902 | STANDARD B | yellow | 15 | 85 | 4 ml |
| BA R-8903 | STANDARD C | orange | 50 | 284 | 4 ml |
| BA R-8904 | STANDARD D | blue | 150 | 851 | 4 ml |
| BA R-8905 | STANDARD E | grey | 500 | 2,840 | 4 ml |
| BA R-8906 | STANDARD F | black | 2,500 | 14,175 | 4 ml |
| BA R-8951 | CONTROL 1 | green | Refer to QC-Report for expected value and acceptable range! | | 4 ml |
| BA R-8952 | CONTROL 2 | red | | | 4 ml |

Conversion: serotonin [ng/ml] x 5.67 = serotonin [nmol/l]

Content: TRIS buffer with non-mercury preservatives, spiked with defined quantity of serotonin

4.3 Additional materials required but not provided in the kit

- Water (deionized, distilled, or ultra-pure)
- Absorbent material (paper towel)
- Reaction tubes, at least 3 ml, Polypropylene/Polystyrol

4.4 Additional equipment required but not provided in the kit

- Calibrated precision pipettes to dispense volumes between 20 – 500 µl
- Microtiter plate washing device (manual, semi-automated or automated)
- ELISA reader capable of reading absorbance at 450 nm and if possible 620 – 650 nm
- Microtiter plate shaker (shaking amplitude 3 mm; approx. 600 rpm)
- Vortex mixer

5. Sample collection and storage

Serum

Collect blood by venipuncture, allow to clot, and separate serum by centrifugation according to manufacturer's instructions at room temperature. Do not centrifuge before complete clotting has occurred. Patients receiving anticoagulant therapy may require increased clotting time. Serum serotonin levels may fluctuate throughout the day. Therefore, the blood sample should always be taken at the same time of day. Traumatic vascular access can drastically increase serotonin levels.

When in doubt, it is recommended that hemolytic, icteric, and lipemic samples not be used in the assay (see 2.2.1).

Storage: up to 1 day at 18 – 25 °C; up to 3 days at 2 – 8 °C; storage for a longer period (up to 6 months) at -20 °C.

Repeated freezing and thawing should be avoided.

Always store samples protected from light.

Urine

24-hour urine samples as well as spontaneous urine (second morning urine) can be used for analysis. 24-hour urine: Over a defined period of 24 hours, all urine is collected in a bottle with acid (10 – 15 ml 100% acetic acid) provided for stabilization and the total volume is noted for evaluation of the results. During the collection period, the collected sample must always be stored in a cool place protected from light (2 – 8 °C).

Spontaneous urine (second morning urine): stabilized with 10 µl 100% acetic acid per 1 ml of urine sample can be used. Always store samples protected from light. A creatinine determination for normalization is required.

When stabilizing urine, consider the acidity (see 2.2.1).

Storage: up to 1 day at 18 – 25 °C; up to 3 days at 2 – 8 °C; storage for a longer period (up to 6 months) at -20 °C.

Repeated freezing and thawing should be avoided.

6. Test procedure

Allow all reagents and samples to reach room temperature and mix thoroughly by gentle inversion before use. Number the reaction tubes and microwell plates (microtiter strips which are removed from the frame for usage should be marked accordingly to avoid any mix-up). Duplicate determinations are recommended.

The binding of the antisera and of the enzyme conjugate and the activity of the enzyme are temperature dependent. The higher the temperature, the higher the absorption values will be. Varying incubation times will have similar influences on the absorbance. The optimal temperature during the enzyme immunoassay is between 20 – 25 °C.

⚠ The use of a microtiter plate shaker with the following specifications is mandatory: shaking amplitude 3 mm; approx. 600 rpm. Shaking with differing settings might influence the results.

⚠ Do not exceed the temperature during the enzyme immunoassay of 20 – 25 °C and the prescribed incubation times. Too high temperature during the enzyme immunoassay and too long incubation times might influence the results.

⚠ Carry out the washing steps thoroughly! Poor washing might influence the results.

6.1 Preparation of reagents and further notes

Wash Buffer

Dilute the 20 ml Wash Buffer Concentrate **WASH-CONC 50x** with water to a final volume of 1000 ml.
Storage: 2 months at 2 – 8 °C

Serotonin Microtiter Strips

In rare cases residues of the blocking and stabilizing reagent can be seen in the wells as small, white dots or lines. These residues do not influence the quality of the product.

Acylation Reagent

The Acylation Reagent (BA E-6612) has a freezing point of 18.5 °C. To ensure that the Acylation Reagent is liquid when being used, it must be ensured that the Acylation Reagent has reached room temperature and forms a homogeneous, crystal-free solution before being used.

If more than 3 ml are needed, pool the contents of the individual vials **ACYL-REAG** and mix thoroughly.

6.2 Preparation of samples – Acylation

- | |
|---|
| 1. Pipette 20 µl of the standards, controls, and samples into the respective reaction tubes. |
| 2. Add 500 µl ACYL-BUFF to all tubes. |
| 3. Add 50 µl of ACYL-REAG to all tubes. |
| 4. Mix the reaction tubes thoroughly (vortex) and incubate for 15 min at RT (20 – 25 °C). |
| 5. Add 500 µl water to all tubes and mix thoroughly (vortex). |

⚠ Take 20 µl of the acylated standards, controls, and samples for the Serotonin ELISA.

6.3 Serotonin ELISA

1. Pipette **20 µl** of the **acylated standards, controls, and samples** into the appropriate wells of the **W SER 5-HIAA**.
 2. Pipette **50 µl** of the **SER-AS** into all wells.
 3. Incubate **60 min at RT** (20 – 25 °C) on a **shaker** (approx. 600 rpm).
 4. Discard or aspirate the content of the wells. Wash the plate **4 times** by adding **300 µl** of **Wash Buffer, discarding** the content and **blotting dry each time** by tapping the inverted plate on absorbent material.
 5. Pipette **100 µl** of the **CONJUGATE** into each well.
 6. Incubate **for 30 min at RT** (20 – 25 °C) on a **shaker** (approx. 600 rpm).
 7. Discard or aspirate the content of the wells. Wash the plate **4 times** by adding **300 µl** of **Wash Buffer, discarding** the content and **blotting dry each time** by tapping the inverted plate on absorbent material.
 8. Pipette **100 µl** of the **SUBSTRATE** into each well.
 9. Incubate **for 25 ± 5 min at RT** (20 – 25 °C) on a **shaker** (approx. 600 rpm).
- ⚠ Avoid exposure to direct sunlight!**
10. Add **100 µl** of the **STOP-SOLN** to all wells and shake the microtiter plate shortly.
 11. **Read** the absorbance of the solution in the wells within 10 min, using a microtiter plate reader set to **450 nm** (if available a reference wavelength between 620 nm and 650 nm is recommended).

7. Calculation of results

| Measuring range | Serotonin | |
|-----------------|-----------|-----------------|
| | Serum | 8 – 2,170 ng/ml |
| | Urine | 8 – 2,027 ng/ml |

The standard curve is obtained by plotting the absorbance readings (calculate the mean absorbance) of the standards (linear, y-axis) against the corresponding standard concentrations (logarithmic, x-axis) using a concentration of 0.001 ng/ml for Standard A (this alignment is mandatory because of the logarithmic presentation of the data).

Use non-linear regression for curve fitting (e. g., 4-parameter, marquardt).

⚠ This assay is a competitive assay. This means: the OD-values are decreasing with increasing concentrations of the analyte. OD-values found below the standard curve correspond to high concentrations of the analyte in the sample and must be reported as being positive.

The concentrations of the samples and controls can be read directly from the standard curve. Samples found with concentrations higher than the highest standard (Standard F) should be diluted accordingly with Standard A and must be re-assayed.

The total amount of **Serotonin** excreted in urine during 24h is calculated as following:

$$\mu\text{g}/24\text{h} = \mu\text{g/l} \times \text{l}/24\text{h}$$

The amount of **Serotonin** normalized to creatinine is calculated as following:

$$\mu\text{g/g creatinine} = \text{ng/ml (serotonin)} / \text{creatinine (mg/dl)} \times 100$$

Conversion

$$\text{Serotonin [ng/ml]} \times 5.67 = \text{serotonin [nmol/l]}$$

7.1 Expected reference value


It is strongly recommended that each laboratory should determine its own reference values. As a basis for the internal reference range determination, the following number of samples for the respective parameters was considered: 24-hour urine n = 194, spontaneous urine (second morning urine) n = 81, serum n = 80. Expected reference ranges were determined in an internal study by testing samples from an apparently healthy European population (95% reference interval).

| | Serotonin |
|--|--|
| Reference range 24-hour urine | 9 – 193 µg/24h 24 – 124 µg/g creatinine |
| Reference range Spontaneous urine (Second morning urine) | 30 – 129 µg/g creatinine |
| Reference range serum | 20 – 206 ng/ml |

Values significantly outside the reference range should be assessed by a doctor.

7.2 Typical standard curve

 Example do not use for calculation!

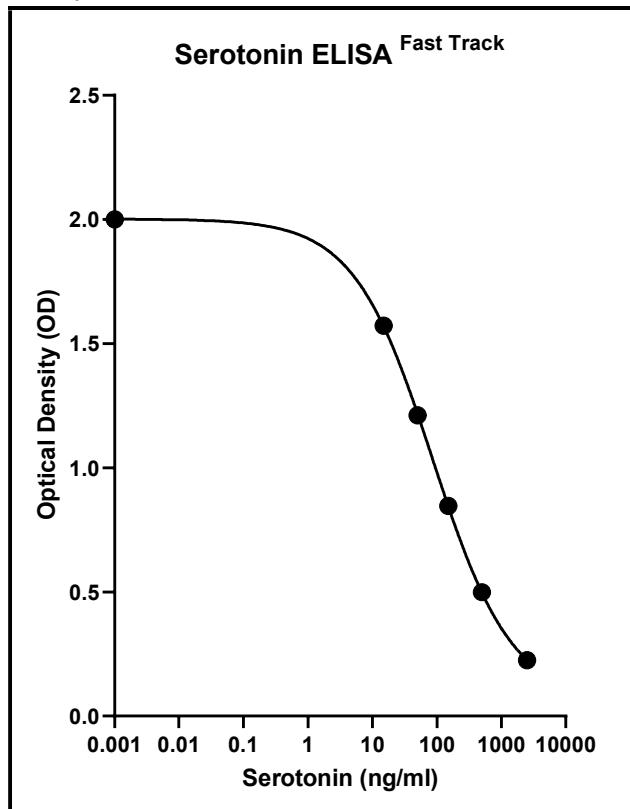8. Quality control

It is recommended to use control samples according to national regulations. Use controls at both normal and pathological levels. Commercially obtained control samples should be treated like unknown samples. Control samples should fall within established confidence limits. The confidence limits of the kit controls are printed on the QC-Report.

9. Assay characteristics

9.1 Performance data

Analytical Sensitivity was determined according to the CLSI Standard EP17-A2 Vol. 32 No. 8. For the determination of the analytical sensitivity, 5 blank samples and 5 low level samples in 2 kit lots in 4 replicates per sample were determined. This resulted in 60 results blank and 60 results low level per lot.

| Analytical Sensitivity | Serotonin |
|-------------------------------|------------------|
| Limit of Blank (LOB) | 2.9 ng/ml |
| Limit of Detection (LOD) | 5.9 ng/ml |
| Limit of Quantification (LOQ) | 8.0 ng/ml |

| Analytical Specificity (Cross Reactivity) | |
|--|-----------------------------|
| Substance | Cross Reactivity (%) |
| Tryptamine | 0.171 |
| Melatonin | < 0.1 |
| 5-Hydroxyindole acetic acid | < 0.1 |
| Phenylalanine | < 0.1 |
| Histidine | < 0.1 |
| Tyramine | < 0.1 |
| 5-Hydroxytryptophan | < 0.1 |

The method comparison was conducted according to the CLSI standard CLSI EP09c 3rd ed.

| | | |
|--|-------|---|
| Method comparison ELISA vs. XLC-MS/MS | Serum | $y = 0.99x - 9.2; r^2 = 0.996; n = 100$ |
| | Urine | $y = 0.9x - 20.7; r^2 = 0.988; n = 97$ |

The precision of the intra- and inter-assay variation was investigated by determining the concentration of 6 serum samples and 6 urine samples in two runs per day in each 2 replicates over 20 days (according to the CLSI standard EP05-A3 Vol. 34 No.13).

| Precision | | | | | |
|--------------------|-----------------------|--------|--------------------|-----------------------|--------|
| Intra-Assay | | | Inter-Assay | | |
| Serum | | | Serum | | |
| Sample | Mean \pm SD [ng/ml] | CV [%] | Sample | Mean \pm SD [ng/ml] | CV [%] |
| 1 | 11.8 \pm 2.1 | 17.6 | 1 | 11.8 \pm 3.3 | 28.2 |
| 2 | 61.6 \pm 5.2 | 8.4 | 2 | 61.6 \pm 7.7 | 12.5 |
| 3 | 102 \pm 8.6 | 8.5 | 3 | 102 \pm 12.3 | 12.1 |
| 4 | 227 \pm 15.5 | 6.8 | 4 | 227 \pm 23.0 | 10.1 |
| 5 | 493 \pm 25.2 | 5.1 | 5 | 493 \pm 55.7 | 11.3 |
| 6 | 1,792 \pm 109 | 6.1 | 6 | 1,792 \pm 165 | 9.2 |
| Urine | | | Urine | | |
| Sample | Mean \pm SD [ng/ml] | CV [%] | Sample | Mean \pm SD [ng/ml] | CV [%] |
| 1 | 18.1 \pm 2.0 | 11.3 | 1 | 18.1 \pm 4.0 | 22.2 |
| 2 | 55.2 \pm 4.0 | 7.3 | 2 | 55.2 \pm 6.4 | 11.7 |
| 3 | 153 \pm 9.1 | 5.9 | 3 | 153 \pm 14.6 | 9.5 |
| 4 | 240 \pm 11.4 | 4.8 | 4 | 240 \pm 21.9 | 9.1 |
| 5 | 498 \pm 29.3 | 5.9 | 5 | 498 \pm 44.5 | 8.9 |
| 6 | 1,798 \pm 120 | 6.7 | 6 | 1,798 \pm 221 | 12.3 |

Recovery was determined according to the CLSI standard EP 34 1st ed.

| Recovery | | | |
|-------------------------------------|-----------------------|----------|-----------|
| | Range [ng/ml] | Mean [%] | Range [%] |
| Serum | 49.4 – 1,046 | 98 | 84 – 112 |
| Urine | 10.0 – 1,023 | 91 | 82 – 98 |
| Linearity of sample dilution | | | |
| | Serial dilution up to | Mean [%] | Range [%] |
| Serum | 1:64 | 103 | 93 – 113 |
| Urine | 1:64 | 98 | 88 – 111 |

The linearity within the measuring range was determined according to the CLSI standard CLSI EP06-Ed2. The linearity is given if the determined value does not deviate by more than 20% from the forecast value.

| Linear range | |
|---------------------|------------------|
| Serum | 18 – 2,170 ng/ml |
| Urine | 20 – 2,027 ng/ml |

| Lot-to-Lot | | | |
|----------------------------|--------|-----------------------|--------|
| | Sample | Mean \pm SD [ng/ml] | CV [%] |
| Serotonin in urine (n = 3) | 1 | 118 \pm 8.5 | 7.2 |
| | 2 | 549 \pm 47.3 | 8.6 |
| Serotonin in serum (n = 3) | 1 | 90.3 \pm 2.0 | 2.2 |
| | 2 | 755 \pm 52.2 | 6.9 |

9.2 Metrological Traceability

The values assigned to the standards and controls of the Serotonin ELISA ^{Fast Track} are traceable to SI Units by calibrated weighing with quality-controlled analyte.

| Standards und Controls | Uncertainty [%] |
|-------------------------------|-----------------|
| | 1.2% |

| Serotonin ELISA ^{Fast Track} | | |
|--|-----------------------|---------------------------------|
| | Concentration [ng/ml] | Expanded Uncertainty [%] k = 2* |
| Serum | 61.6 | 25.1 |
| | 227 | 20.3 |
| Urine | Concentration [ng/ml] | Expanded Uncertainty [%] k = 2* |
| | 18.1 | 44.5 |
| | 55.2 | 23.5 |
| | 153 | 19.2 |
| | 240 | 18.4 |
| | 498 | 18.0 |
| | 1,798 | 24.7 |

* This defines an interval about the measured result that will include the true value with a probability of 95%.

10. References/Literature

1. Bieger, W.P., NeuroStress Guide. 2011.
2. Chojnacki, C., et al., Evaluation of serotonin and dopamine secretion and metabolism in patients with irritable bowel syndrome. *Pol Arch Intern Med*, 2018. 128(11): p. 711 – 713.
3. Huang, H., Z. Chen, and X. Yan, Simultaneous determination of serotonin and creatinine in urine by combining two ultrasound-assisted emulsification microextractions with on-column stacking in capillary electrophoresis. *J Sep Sci*, 2012. 35(3): p. 436 – 44.
4. Piešťanský, J., K. Maráková, and P. Mikuš, Two-Dimensional Capillary Electrophoresis with On-Line Sample Preparation and Cyclodextrin Separation Environment for Direct Determination of Serotonin in Human Urine. *Molecules*, 2017. 22(10).
5. Lindström, M., et al., Comparison of serum serotonin and serum 5-HIAA LC-MS/MS assays in the diagnosis of serotonin producing neuroendocrine neoplasms: A pilot study. *Clin Chim Acta*, 2018. 482: p. 78 – 83.
6. Ren, C., et al., Low levels of serum serotonin and amino acids identified in migraine patients. *Biochem Biophys Res Commun*, 2018. 496(2): p. 267 – 273.
7. Moriarty, M., et al., Development of an LC-MS/MS method for the analysis of serotonin and related compounds in urine and the identification of a potential biomarker for attention deficit hyperactivity/hyperkinetic disorder. *Anal Bioanal Chem*, 2011. 401(8): p. 2481 – 93.
8. Nichkova, M.I., et al., Evaluation of a novel ELISA for serotonin: urinary serotonin as a potential biomarker for depression. *Anal Bioanal Chem*, 2012. 402(4): p. 1593 – 600.
9. Holck, A., et al., Plasma serotonin levels are associated with antidepressant response to SSRIs. *J Affect Disord*, 2019. 250: p. 65 – 70.
10. Jaworek, A.K., et al., Depression and Serum Content of Serotonin in Adult Patients with Atopic Dermatitis. *Adv Exp Med Biol*, 2020. 1271: p. 83 – 88.
11. Shu, B., et al., Serotonin and YAP/VGLL4 Balance Correlated with Progression and Poor Prognosis of Hepatocellular Carcinoma. *Sci Rep*, 2018. 8(1): p. 9739.
12. Mahato, K., et al., Novel electrochemical biosensor for serotonin detection based on gold nanorattles decorated reduced graphene oxide in biological fluids and in vitro model. *Biosens Bioelectron*, 2019. 142: p. 111502.
13. Yu, R., High Serum Serotonin Test Results Caused by Traumatic Vascular Access Due to Difficult Veins. *Pancreas*, 2019. 48(6): p. e51 – e53.
14. Mamdouh, F., et al., Serum Serotonin as a Potential Diagnostic Marker for Hepatocellular Carcinoma. *J Interferon Cytokine Res*, 2019. 39(12): p. 780 – 785.

For updated literature or any other information please contact your local supplier.

11. Changes

| Version | Release Date | Chapter | Change |
|---------|--------------|---------|--|
| 17.1 | 2022-06-14 | All | Assay Revision due to IVDR regulation (EU) 2017/746. Due to this revision all chapters were revised and changed. |
| 18.0 | 2023-03-20 | 7.1 | Sentence: "Values significantly outside the reference range should be assessed by a doctor." added |

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10. Referenzen/Literatur

1. Bieger, W.P., NeuroStress Guide. 2011.
2. Chojnacki, C., et al., Evaluation of serotonin and dopamine secretion and metabolism in patients with irritable bowel syndrome. *Pol Arch Intern Med*, 2018. 128(11): p. 711 – 713.
3. Huang, H., Z. Chen, and X. Yan, Simultaneous determination of serotonin and creatinine in urine by combining two ultrasound-assisted emulsification microextractions with on-column stacking in capillary electrophoresis. *J Sep Sci*, 2012. 35(3): p. 436 – 44.
4. Piešťanský, J., K. Maráková, and P. Mikuš, Two-Dimensional Capillary Electrophoresis with On-Line Sample Preparation and Cyclodextrin Separation Environment for Direct Determination of Serotonin in Human Urine. *Molecules*, 2017. 22(10).
5. Lindström, M., et al., Comparison of serum serotonin and serum 5-HIAA LC-MS/MS assays in the diagnosis of serotonin producing neuroendocrine neoplasms: A pilot study. *Clin Chim Acta*, 2018. 482: p. 78 – 83.
6. Ren, C., et al., Low levels of serum serotonin and amino acids identified in migraine patients. *Biochem Biophys Res Commun*, 2018. 496(2): p. 267 – 273.
7. Moriarty, M., et al., Development of an LC-MS/MS method for the analysis of serotonin and related compounds in urine and the identification of a potential biomarker for attention deficit hyperactivity/hyperkinetic disorder. *Anal Bioanal Chem*, 2011. 401(8): p. 2481 – 93.
8. Nichkova, M.I., et al., Evaluation of a novel ELISA for serotonin: urinary serotonin as a potential biomarker for depression. *Anal Bioanal Chem*, 2012. 402(4): p. 1593 – 600.
9. Holck, A., et al., Plasma serotonin levels are associated with antidepressant response to SSRIs. *J Affect Disord*, 2019. 250: p. 65 – 70.
10. Jaworek, A.K., et al., Depression and Serum Content of Serotonin in Adult Patients with Atopic Dermatitis. *Adv Exp Med Biol*, 2020. 1271: p. 83 – 88.
11. Shu, B., et al., Serotonin and YAP/VGLL4 Balance Correlated with Progression and Poor Prognosis of Hepatocellular Carcinoma. *Sci Rep*, 2018. 8(1): p. 9739.
12. Mahato, K., et al., Novel electrochemical biosensor for serotonin detection based on gold nanorattles decorated reduced graphene oxide in biological fluids and in vitro model. *Biosens Bioelectron*, 2019. 142: p. 111502.
13. Yu, R., High Serum Serotonin Test Results Caused by Traumatic Vascular Access Due to Difficult Veins. *Pancreas*, 2019. 48(6): p. e51-e53.
14. Mamdouh, F., et al., Serum Serotonin as a Potential Diagnostic Marker for Hepatocellular Carcinoma. *J Interferon Cytokine Res*, 2019. 39(12): p. 780 – 785.

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SYMBOLS USED WITH DEMEDITEC ASSAYS

| Symbol | English | Deutsch | Française | Espanol | Italiano |
|-----------------------|------------------------------------|--|--|---|-------------------------------------|
| | European Conformity | CE-Konformitäts-kennzeichnung | Conforme aux normes européennes | Conformidad europea | Conformità europea |
| | Consult instructions for use | Gebrauchsanweisung beachten | Consulter les instructions d'utilisation | Consulte las Instrucciones | Consultare le istruzioni per l'uso |
| | In vitro diagnostic device | In-vitro-Diagnostikum | utilisation Diagnostic in vitro | Diagnóstico in vitro | Per uso Diagnostica in vitro |
| | For research use only | Nur für Forschungszwecke | Seulement dans le cadre de recherches | Sólo para uso en investigación | Solo a scopo di ricerca |
| | Catalogue number | Katalog-Nr. | Référence | Número de catálogo | No. di catalogo |
| | Lot. No. / Batch code | Chargen-Nr. | No. de lot | Número de lote | Lotto no |
| | Contains sufficient for <n> tests/ | Ausreichend für "n" Ansätze | Contenu suffisant pour "n" tests | Contenido suficiente para <n> ensayos | Contenuto sufficiente per "n" saggi |
| | Note warnings and precautions | Warnhinweise und Vorsichtsmaßnahmen beachten | Avertissements et mesures de précaution font attention | Tiene en cuenta advertencias y precauciones | Annoti avvisi e le precauzioni |
| | Storage Temperature | Lagerungstemperatur | Température de conservation | Temperatura de conservacion | Temperatura di conservazione |
| | Expiration Date | Mindesthaltbarkeits-datum | Date limite d'utilisation | Fecha de caducidad | Data di scadenza |
| | Legal Manufacturer | Hersteller | Fabricant | Fabricante | Fabbricante |
| <i>Distributed by</i> | Distributed by | Vertrieb durch | Distribution par | Distribución por | Distribuzione da parte di |
| | Version | Version | Version | Versión | Versione |
| | Single-use | Einmalverwendung | À usage unique | Uso único | Uso una volta |